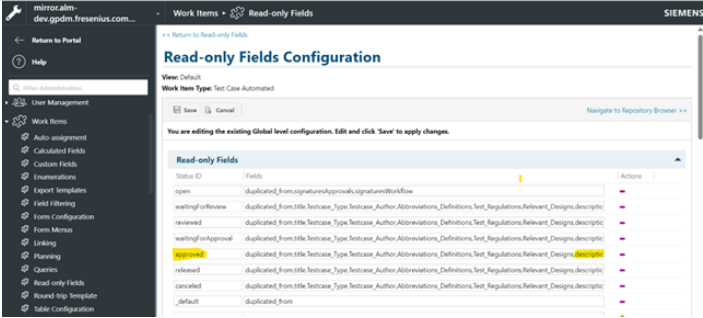
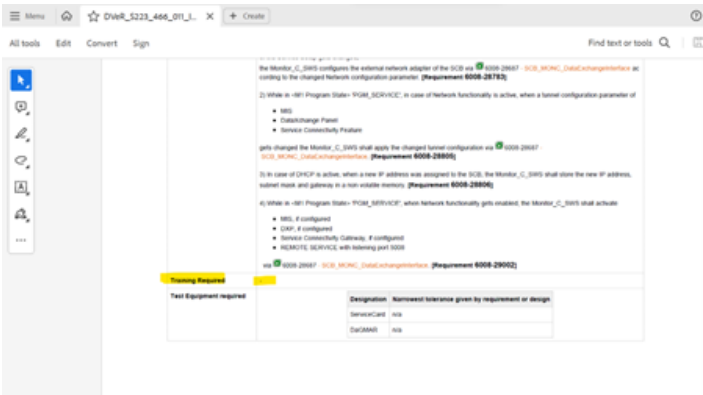
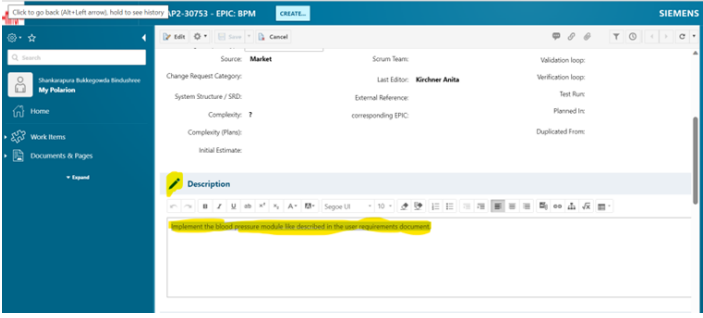


How to handle formatting issues in Document

Step	Action	Screenshot
1	Note: If the work items are in "draft" status, possible to edit the attribute values. If work items are in other status eg: approved, released etc..., not possible to edit the attribute values of the work item directly as in those status attribute values are read-only. Steps to be followed by AMS Team to enable the attribute values. 1. Export the document as PDF. 2. Identify the work item which is causing the formatting issue. 3. Then check the type and status of the work item, if the status is not in draft and not able to edit attribute values. 4. Go to Administration -> Work items -> Read-only fields -> select the type of work item 5. Delete the particular attribute from the particular status of the WI (Eg: description) then save the read-only fields. 6. After correcting the formatting issues in work items, again add the same fields with the particular status of the WI. Save the Read-only fields.	
2	Case 1: After export of a document, some fields does not display content instead it shows "-". Root Cause: In work item, even the "Training Required" field has the content. After exporting the document it shows "-" instead of content User copy pasted the content with some different format that causes the formatting issue. Solution: • Identify the work items which is causing the issue. • It is usually corrected by copy the content to word document then do merge formatting and pasting the content again by editing the required fields (the field which is causing formatting issue).	 

3

Case 2:

After PDF export, font size of the text is too small to be read.



Root Cause:

Ex:1

while observing the document itself able to identify, whether it has formatting issue or not.

In the document, scroll bar should not be there at the bottom. If it has scroll bar, there is a formatting issue. Need to identify which work item is causing this issue.

Identified the WI having issue. It has a very big image; this image needs to be re-sized and remove the extra space then save the work item.

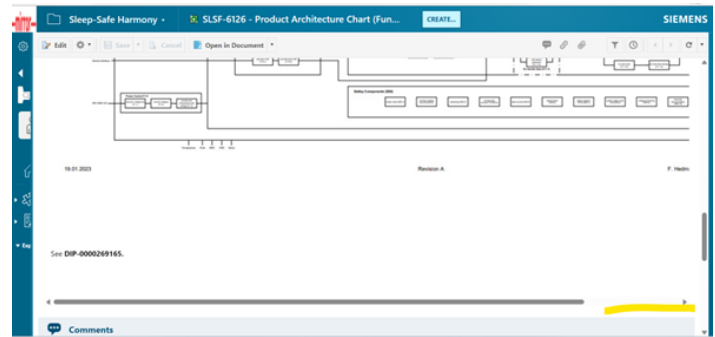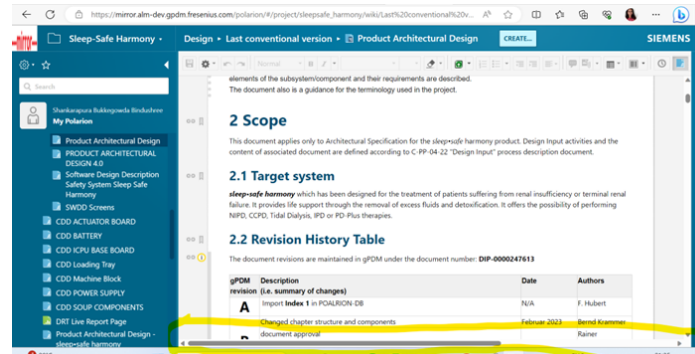

4

EX:2

In the PDF, the work item which exceeding the space causing the formatting issue.

Solution:

- Identify the work items which is causing the issue.
- Edit the description by copy pasting the content, if required Re-size the images and save the work item
- Refresh the document and export the PDF.



5

Case 3:

After PDF Export, some text passages highlighted with yellow colour.

Root Cause:

Formatting issue with the content which user used in the description field.

Solution:

- Identify the work items which is causing the issue.
- Edit the description field, Select the text which is highlighted with yellow color.
- Copy the text and paste in word document and do merge formatting
- Again, copy the content from Word document to description field of the Work item.
- Save the work item.
- Refresh the document and export PDF

- Close V43.
- Stop Air Separation Pump.
- Stop UF pump.
- Set Balancing Chamber to all valves off state.

Depending on the cleaning phases, the following actions will be taken additionally to the general behavior described above:

- Empty rinse chambers and air separation pots (part of rinse similar programs, see 4AP2-65733):
 - Continue flow to next overlap time:
 - Open V42 at 46, 47, 53 and 54 timeslices after the change to overlap time if filling mode is active.
 - Control V43 as if there is no interruption.
 - Instant stop of flow:
 - Close V42.
- Fill rinse chambers and air separation pots (part of rinse similar programs, see 4AP2-65726):
 - Instant stop of flow:
 - Open V42.
 - (State Continue flow to next overlap time not relevant since no balancing chamber switches during this subphase)
- Periodic rinsing of hydraulics (part of rinse similar programs, see 4AP2-65735),
Periodic filling of hydraulics (part of rinse similar programs, see 4AP2-65736),
Degreasing (see 4AP2-35452),
Mandatory rinsing (see 4AP2-35454),
Heat disinfection (see 4AP2-35450):
 - Continue flow to next overlap time:
 - Open V42 instead of the valve to be opened in the mixing circuit (see additional information below),
 - Control V43, air separation pump and ultrafiltration pump as if there is no interruption.
 - Instant stop of flow:
 - Close V42.

6

Case 4:

Ex1:


After exporting the document columns of tables are not aligned properly

Root Cause:

- In WorkItem.xml file and found ghost html object (a container consuming space but is empty) and as a result the table columns division is not even and table appears non-aligned in the exported PDF.
- This was caused by copy paste action performed by user.

Solution:

Removing the highlighted line from code resolved the misalignment in the exported pdf.

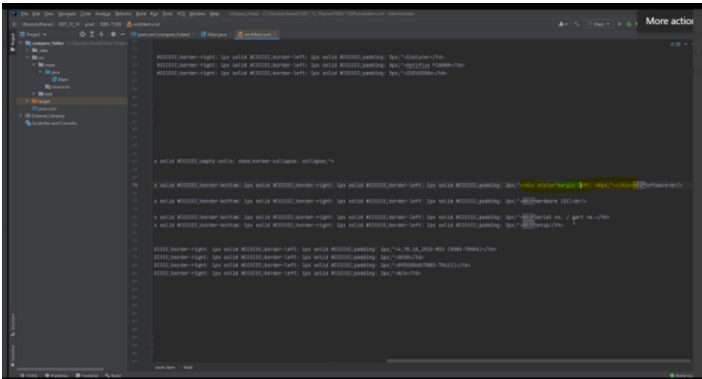

7

Ex2:

After exporting the document columns of tables are not aligned properly.

Root Cause:

Formatting issue with the content which user used in the Test Step field.

Solution:

- Identify the work items which is causing the issue.
- In Test Steps field, delete all the content inside the table and save the work item.
- Copy paste the same content again into the table and click on "Remove formatting of the selected text" button.
- Save the work item.
- Re-run the test run.
- Export of the document

Protocol execution information			
Step ID	Step Description	Expected Result	Actual Result
or internal use only - Page 14 of 29 Revision 3139204			
VIR-0000207990 GDL Physical CO2 Barrier Verification Report Configuration Template: 2.0.12.3			
1	Physical Testing - Primary Bag CO2 Barrier Testing	Acceptance Criteria: • The bag foil shall have sufficient barrier characteristics to minimize water & CO2 loss. • Primary bag foil shall limit CO2 loss	12 Samples were evaluated • The bag foil has sufficient barrier characteristics to minimize water & CO2 loss. • Primary bag foil limits CO2 loss over shelf life • CO2 Leak Limit is < 20cm³/m² dba

Test Step Number	Step Description	Expected Result
1	Physical Testing - Primary Bag CO2 Barrier Testing Sample Size: 10 Requirements: • P295-16932 - Primary bag foil - permeability • P295-17277 - Primary Bag Foil - Vapor Barrier Test Method: Physical • Test each sample using the specified Test Method and Equipment	Acceptance Criteria: • The bag foil shall have sufficient barrier characteristics to minimize water & CO2 loss. • Primary bag foil shall limit CO2 loss • Limit < 20cm³/m² dba

Test Step Number	Step Description	Expected Result
1	Physical Testing - Primary Bag CO2 Barrier Testing Sample Size: 10 Requirements: • P295-16932 - Primary bag foil - permeability • P295-17277 - Primary Bag Foil - Vapor Barrier Test Method: Physical • Test each sample using the specified Test Method and Equipment	Acceptance Criteria: • The bag foil shall have sufficient barrier characteristics to minimize water & CO2 loss. • Primary bag foil shall limit CO2 loss • Limit < 20cm³/m² dba

8

Case 5:

1. AMS Team have received a couple of requests of formatting issue from end users.
2. When user do a PDF export, some fields (Training required) in the document will show "-" instead of content.

the Monitor_C_3045 configures the external network adapter of the SCS via 6209-28687 - SCS_MNC_DataExchangeInterface as	
2) While in -HMI Program State- PGM_SERVICE, in case of network functionality is active, when a tunnel configuration parameter of	
• SCS	
• DataExchange Panel	
• Service Connectivity Protocol	
gets changed the Monitor_C_3045 shall apply the changed tunnel configuration via 6209-28687	
SCS_MNC_DataExchangeInterface (Requirement 6209-28687)	
3) In case of SCS-F is active, when a new IP address was assigned to the SCS, the Monitor_C_3045 shall store the new IP address,	
related mode and gateway in a non-volatile memory. (Requirement 6209-28686)	
4) While in -HMI Program State- PGM_SERVICE, when network functionality gets enabled, the Monitor_C_3045 shall activate	
• SCS, if configured	
• SCS-F, if configured	
• Service Connectivity Gateway, if configured	
• REMOTE_SERVICE, with relevant port SCS	
via 6209-28687 - SCS_MNC_DataExchangeInterface (Requirement 6209-28682)	
Training Required	
Test Equipment required	Designation: Network tolerance given by requirement or design
	NetworkCard - N/A
	Car2045 - N/A

3. AMS Team knows that it is a formatting issue and usually corrected by copy pasting the content by editing the required fields (the field which is causing formatting issue).

4. In specific cases with the formatting issue,

Design Verification Interim Summary or Design Verification report.

User has to re-run the Test run, because in the report the revision of the workitem when it was executed in the test run will be fetched.

5. Steps to be followed:

1. Identify the work items which is causing the formatting issue.

2. Re-test the workitems with the same test run.

3. After re-run, if the user do the export of the document, in the specified fields they will be able to see the content instead of "-".

6. Sometimes, users are not ready to re-execute the Test run. There is a work around to edit the testrun.xml file and update the latest revision of the workitem in the testrun.xml file.

Note: It is not a recommended approach to update the testrun.xml file directly from the AMS Team, however, this can be done on request from the project team.



7. If the user insists on making the changes and does not want to re-run the Test run, then in that case the following steps have to be performed.

- Inform the user to involve respective QA (Quality Assurance person) to the ticket to maintain transparency.
- So, the QA knows what changes are done.

8. AMS Support also double check after re-execute of the Test Run, if the Test Protocol work item still has the same revision/content as at the moment the Test Protocol was executed.

AMS Team should follows the below steps:

- Open the document which has formatting issue.
- Compare the work item's Revision of Test Protocol Approval in the test run with the history of Work item's Revision of Test Protocol Approval.

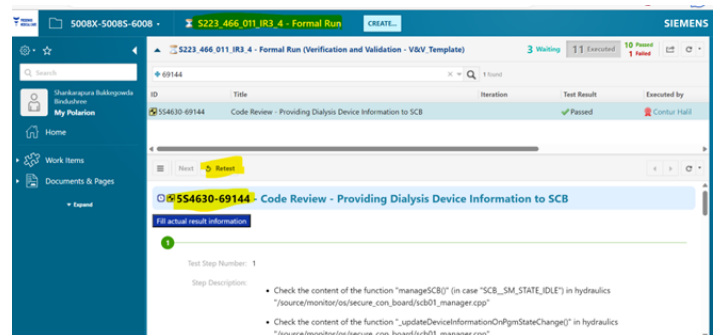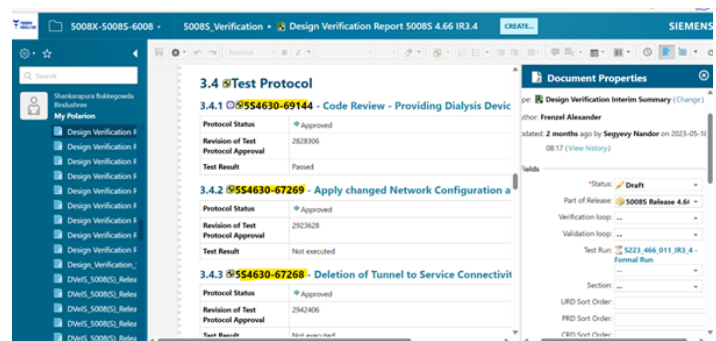

- Both the Revisions should be same, so that checking is completed.

The screenshot displays the Siemens Polarion interface for a project named '5008X-5008S-6008'. The main content area shows three code review entries for item '554630-69143'.

Section	Protocol Status	Revision of Test Protocol Approval	Test Result
3.4.8 554630-69143 - Code Review - Triggering Switch On/Off Instruction Detection System	Approved	3040966	Passed
3.4.9 554630-69143 - Code Review - Providing Current Program State to SCB	Approved	2628306	Passed
3.4.10 554630-67400 - Displaying Error Message, if SCB is not Responsive	Approved	2937879	Passed

On the right side, a sidebar shows document details for 'Design Verification Report 5008S 4.66 IR3.4', including the author 'Frenzel Alexander' and update history.

The screenshot shows the history of item '554630-69143' in the Siemens Polarion interface. The history table tracks updates, status changes, and approvals.

Field	Removed	Added
Updated By: Contur Hall	Date: 2023-06-23 07:42, Revision: 3040998	
Status	Specified	Approved
Last Editor	Hedmann Frank	Contur Hall
Revision of Test Protocol Approval	3040966	
Updated By: Hedmann Frank	Date: 2023-06-23 07:35, Revision: 3040966	
Comments		#6 Item was e-signed by Hedmann Frank on 2023-06-23 07:35. Approval verdict: Approved
Approvals	Hedmann Frank: Waiting	Hedmann Frank: Approved
Last Editor	Contur Hall	Hedmann Frank
signatures Approvals		Item was e-signed by Contur Hall (h5850u) on 2023-06-23 07:24 during approval. Verdict: Approved
		Item was e-signed by Hedmann Frank (h5943) on 2023-06-23 07:35 during approval. Verdict: Approved
Updated By: Contur Hall	Date: 2023-06-23 07:24, Revision: 3040998	